Otto Bock®

QUALITY FOR LIFE

Motus

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Instructions for Use for the Motus Adaptive Wheelchair

Contents	Page
1 General Information	41
1.1 Preface	41
1.2 Intended Use	41
1.3 Field of Application	41
1.4 CE Conformity	42
1.5 Liability	42
1.6 Service	42
2 Safety Instructions	43
2.1 Explanation of Symbols	43
2.2 General Safety Instructions	43
3 Delivery / Preparing the Wheelchair for Use	47
4 Transportation	48
4.1 Transfer	49
5 Options	51
5.1 Anti-tipper	51
5.2 Transport Wheels	51
5.3 Crutch Holder	52
5.4 Tip-assist	52
5.5 Plug-on Side Panels	52
5.6 Desk Side Panels	52
5.7 Plug-on Wheel Lock Lever Extension	53
5.8 Rear Wheel with Drum Brake	53
5.9 Tray	53
5.10 Elevating Footrest	54
5.11 Push Handles, Height-adjustable	54
5.12 Wheelbase Extension	54
5.13 Angle Adjustable Back	54
5.14 Spoke Protector	54
6 Adjustment / Assembly Instructions	55
6.1 Adjusting Seat Height and Seat Inclination	55
6.2 Changing the Wheelbase	56
6.3 Setting the Wheel Camber and Toe-in	56
6.4 Caster Trailing Adjustment	58
6.5 Adjusting the Wheel Lock	58

<i>6.6 Adjusting the Brake Force - Drum Brake.....</i>	<i>59</i>
<i>6.7 Angle Adjustment of the Footplate</i>	<i>59</i>
<i>6.8 Adjusting the Lower Leg Length</i>	<i>60</i>
<i>6.9 Adjusting the Angle of the Elevating Footrest</i>	<i>60</i>
<i>6.10 Adjusting the Height of the Armrest for Desk Side Panel.....</i>	<i>60</i>
<i>6.11 Adjusting the Back Height.....</i>	<i>61</i>
<i>6.12 Installing the Anti-tipper</i>	<i>62</i>
<i>6.13 Adapting the Back Upholstery</i>	<i>63</i>
<i>6.14 Changing the Back Angle.....</i>	<i>63</i>
<i>6.15 Push Handles, Height-adjustable.....</i>	<i>64</i>
<i>6.16 Attaching the Safety Belt</i>	<i>64</i>
7 Maintenance, Cleaning, and Disinfection	65
<i>7.1 Maintenance</i>	<i>65</i>
<i>7.2 Cleaning and Disinfection</i>	<i>66</i>
8 Changing Tires	66
9 Technical Data (in cm)	68
10 Information on Re-use	70

1 General Information

1.1 Preface

With the purchase of the Motus adaptive wheelchair, you have selected a quality product for versatile use during everyday life, both indoors and outdoors. Before using your wheelchair, it is essential to read the sections “Delivery / Preparing the Wheelchair for Use” and “Safety Instructions”.

The chapter “Options” presents various additional components of the Motus, which can expand its field of application and improve comfort for the wheelchair user. The chapter “Adjustment / Assembly Information” provides an overview of the ways the wheelchair can be customized to your requirements. The model as described in these Instructions for Use is subject to technical changes without notice.

1.2 Intended Use

The Motus adaptive wheelchair is designed solely for individual indoor and outdoor use by persons who are unable to walk or who have a walking impediment, and can be operated by the patient or by another person.

The Motus adaptive wheelchair may only be combined with the options mentioned in these Instructions for Use and vice versa. Otto Bock assumes no liability for combinations with medical products and/or accessories from other manufacturers outside of the modular system.

1.3 Field of Application

The versatility of this power wheelchair and its modular design make it suitable for patients who have walking impediments or walking disabilities due to:

- Palsies / paralysis (paraplegia, tetraplegia)
- Loss of limbs (lower limb amputation)
- Defective or deformed limbs
- Joint contractures or defects
- Muscle and nerve disorders
- Diseases such as cardiac or circulatory insufficiency, balance disorders or cachexia as well as for geriatric patients who still have usable residual strength in the upper limbs.

The Motus was specially designed for individuals who are able to independently move in a wheelchair.

Fitting considerations:

- Body size and body weight (max. load capacity 125 or 140 kg in connection with double crossbrace)
- Physical and psychological limitations
- Age of the user
- Home and living conditions, and
- Environment.

1.4 CE Conformity

The product meets the requirements of the European Directive for Medical Products 93/42/EEC. Based on the classification criteria for medical products according to Appendix IX of the Directive, the product was classified as Class I. The declaration of conformity was therefore created by Otto Bock with sole responsibility according to appendix VII of the guidelines.

1.5 Liability

The manufacturer's warranty applies only if the product has been used under the conditions and for the purposes described. The manufacturer recommends that the product be used and maintained according to the instructions for use.

The manufacturer is not responsible for damages caused by components and spare parts not approved by the manufacturer. Repairs must be carried out exclusively by authorized dealers or by the manufacturer.

1.6 Service

Service and repairs on the product may only be carried out by authorized dealers. Should any problems arise, please contact the supplier of your product. Any necessary repairs will be made exclusively with authentic Otto Bock spare parts there. Your product requires regular maintenance / regular service (see Chapter 7).

Your authorized Otto Bock dealer:

2 Safety Instructions

2.1 Explanation of Symbols

 WARNING	Warnings regarding possible risks of severe accident or injury.
 CAUTION	Warnings regarding possible risks of accident or injury.
 NOTICE	Warnings regarding possible technical damage.
 INFORMATION	Information regarding operation. Information for service personnel.

2.2 General Safety Instructions

INFORMATION

Please read the Instructions for Use first! Before using the product, you should become familiar with the handling, function and use of the product. You may compromise the safety of your child if you do not observe these instructions.

However, it is impossible to protect against and avoid all unforeseeable circumstances and situations.

WARNING

Use as a seat for transportation in a vehicle for the disabled

Risk of injury as a result of unauthorized use in a vehicle for the disabled. We recommend that, wherever and whenever possible, Motus users transfer to the seats installed in the motor vehicle and use the corresponding vehicle restraint systems, because this is the only way to ensure optimal protection of the passengers in case of an accident. Only this can ensure optimal protection of the passengers in case of an accident. Your Motus adaptive wheelchair is permissible for transporting passengers in motor vehicles when using Otto Bock safety components as well as appropriate restraint systems.

For more information, please refer to our instructions for use manual "Using your wheelchair/ mobility base with seating shell or buggy as a seat for transportation in motor vehicles", order no. 646D158.

Danger when overcoming stairs/obstacles

WARNING

Risk of falling in case of lacking assistance. Do not ascend or descend stairs without the assistance of attendants. Make use of ramps or elevators whenever available. If wheelchair-friendly access is missing, two attendants must carry the wheelchair over the obstacle.

WARNING

Risk of falling as a result of incorrectly set anti-tipper. If only one attendant is available when ascending or descending stairs, an incorrectly set anti-tipper (if mounted at all) can lead to severe falls. Adjust the anti-tipper so that it does not come in the way of the steps during transport. Afterwards, swing the anti-tipper back to its operational position.

WARNING

Risk of falling as a result of lifting the wheelchair incorrectly. Attendants must lift the wheelchair only by parts that are firmly attached (e.g. not at the footrest or rear wheels). If your wheelchair is equipped with height-adjustable push handles, ensure that the clamping levers are firmly tightened.

WARNING

Risk of tipping when using lifting platforms. Ensure that the anti-tipper (if mounted) is outside the danger-area.

Danger when getting into the wheelchair

⚠ CAUTION

Risk of injury as a result of using the footplate for getting into the wheelchair. If possible, the footplates should be flipped up before getting into or out of the wheelchair.

Danger while driving

⚠ CAUTION

Risk of accident when driving without experience. This can lead to falls and other dangerous situations. Therefore, familiarize yourself with the new wheelchair by practicing on even, straightforward terrain first. Get to know how the wheelchair reacts when the center of gravity shifts; for example on slopes or inclines or when clearing obstacles such as steps and curbs. This should be done only with assistance from another person. Using the anti-tipper is strongly recommended for untrained wheelchair users.

⚠ CAUTION

Risk of tipping when ascending too steep slopes. When ascending slopes or ramps and when crossing obstacles on upward slopes, always lean your upper body far forward.

⚠ CAUTION

Risk of tilting or tipping over as a result of centre of gravity shifts. When reaching for objects in front, to the side or behind the wheelchair, do not lean out of the wheelchair too far.

⚠ CAUTION

Risk of tipping when driving against obstacles (steps, curbs) **without braking.** Use your wheelchair properly. Do not „jump“ the chair down from higher surfaces.

⚠ CAUTION

Risk of accident as a result of incorrect use of the wheel lock. All brakes acting on the tires do not serve as operating brake but are only designed as parking brake. The wheel locks must not be used as driving brakes for slowing down the wheelchair, as in extreme cases, the abrupt stopping of the wheelchair can lead to falls.

⚠ CAUTION

Risk of accident if the wheelchair starts rolling. Engage the wheel locks to prevent the wheelchair from moving on uneven ground or during transfers (e.g. into a car).

⚠ CAUTION

Risk of accident as a result of driving in the dark. In the dark, wear light clothing or clothing with reflectors in order to improve your visibility. Ensure that the reflectors installed on the sides and rear of the wheelchair are easily visible. We also recommend installing lights.

⚠ CAUTION

Risk of tipping when descending slopes. Especially on slopes, do not drive without braking and reduce your speed.

Danger as a result of incorrect assembly or settings

⚠ CAUTION

Risk of tipping as a result of incorrectly set wheels. Extreme settings (e.g. rear wheels mounted in the foremost position) combined with an unfavorable body posture can cause the wheelchair to tip even on level ground. Avoid extreme settings. For transfemoral amputees, the rear wheels must be shifted backwards. Use of an anti-tipper is required.

⚠ CAUTION

Risk of tipping as a result of missing or incorrectly mounted anti-tipper. Under no condition should the anti-tipper be used as transport wheels, e.g. when wanting to push the wheelchair while the rear wheels are removed. The anti-tipper must audibly lock in place, before it is able to bear loads. Firm seating must be verified by the user or by an attendant.

⚠ CAUTION

Risk of tipping as a result of changed diameters/mounting positions of the wheels. If the size and position of the casters or the rear wheel size are modified, the casters may flutter when driving at higher speeds. This can lead to blocking of the casters and cause the wheelchair to tip over. For this reason, you should always maintain the factory settings or, if modifications are required, ensure that the wheelchair frame is aligned horizontally (see "Adjustment / Assembly Instructions").

⚠ CAUTION

Risk of falling as a result of unintentional loosening of the rear wheels. Always make sure that the quick-release axles are correctly set on the rear wheel. The rear wheel must not be removable without the button of the quick-release axle being depressed.

Danger as a result of improper use of the tires

⚠ CAUTION

Accident risk due to insufficient or excessive tire pressure! The effectiveness of the wheel lock and the overall driving quality are dependent on adequate air pressure. The wheel locks are only effective with sufficient air pressure and proper adjustment (approx. **5 mm** distance, subject to technical changes). Before starting to use your mobility base with seating shell, check that the tires are inflated correctly. The required air pressure is printed on the side of the tire. For rear wheels, it should be **at least 7.5 bar (350 kPa)**. Properly inflated rear wheel tires and equal air pressure in both tires considerably improve the maneuverability of your wheelchair.

⚠ CAUTION

Risk of accident as a result of bad tires. Insufficient tread depth of the tires reduces adhesion. When using your wheelchair in public traffic, you are subjected to public traffic regulations.

Risk of skin damage

⚠ CAUTION

Possible complications in case of unhealthy skin. Using the wheelchair with unhealthy skin can cause medically indicated complications such as reddening of the skin and pressure sores. Over several hours of use, the skin on the buttocks, back and rear thighs is subjected to considerable strain. Therefore, always check these areas carefully before using the wheelchair.
No liability will be assumed for health injuries caused by use of this wheelchair with damaged skin.

Danger caused by fire/heat and cold

⚠ CAUTION

Risk of burns when near to fire. The seat and back upholstery may catch fire. Keep the wheelchair away from all sources of ignition, especially lit cigarettes.

⚠ CAUTION

Be careful in case of extreme temperatures. The wheelchair can heat up significantly due to solar radiation or in the sauna. In extreme cold, there is a risk of hypothermia.

Risk of hand injuries

⚠ CAUTION

Risk of hand injuries. To avoid hand injuries, do not grasp between the rear wheel and the wheel lock when driving your wheelchair.

⚠ CAUTION

Risk of pinching between wheelchair components. The wheel lock lever and the side or frame sections can represent a risk of crushing with various adjustment options.

⚠ CAUTION

Risk of friction burns. Slowing down from high speeds or when descending longer slopes tends to heat up the hands and fingers, especially if using aluminum push rings. When using your wheelchair outdoors, wear leather gloves. Gloves provide you with a better grip and protect your fingers from dirt and hot metal.

Warnings regarding wheelchair damage

NOTICE

Damage caused by overload. The maximum load capacity of the Motus adaptive wheelchair is 125 kg or 140 kg in connection with double crossbrace.

NOTICE

Damage caused by wear and tear. If the upholstery gets worn or damaged, replace it immediately.

NOTICE

Damage caused by incorrect packaging. Only use original packaging for shipment of the product.

3 Delivery / Preparing the Wheelchair for Use

⚠ CAUTION

Risk of pinching. Only grip the components shown when folding or collapsing the wheelchair.

Your wheelchair is usually delivered fully assembled with the backrest folded and with the rear wheels removed.

A few simple steps are sufficient to prepare it for use:

1. Attach rear wheels. To attach the rear wheel, press the button on the quick-release axle, slide four fingers between the spokes and press the button of the quick-release axle with your thumb (fig. 1). Then attach the rear wheels to the fittings and the axles will automatically engage.
2. While standing next to the wheelchair, tip it slightly toward you. Press apart the crossbrace with flat of the hand (fig. 2). Make sure that the crossbrace has locked in the support brackets (fig. 3).
3. Pull the bib of the back upholstery upwards and attach the bib to the wheelchair's seat upholstery (fig. 4).
4. Insert the seat cushion (fig. 5).
5. Finally, fold down the footplates (fig. 6).

⚠ CAUTION Accident risk. When remounting, be sure the quick-release axle securely locks into the quick-release axle housing! The rear wheel must not be removable when the button is not being pressed!

⚠ CAUTION Risk of tipping. The rear wheels have been pre-positioned according to your indications. Before using the wheelchair, this pre-adjusted position must be checked by the wheelchair user for function and stability against tipping with the assistance of skilled staff.

4 Transportation

Pull the footrest up and remove the seat cushion.

Grab the seat upholstery at the front and the back in the middle and pull it upwards: The wheelchair will fold (fig. 7).

For more manageable transportation, e.g. in cars, the footrests can be swung outwards or be entirely removed in just a few easy steps.

To move or remove the footrest, pull the swivel handle to the rear (fig. 8). The footrest with footplate can then be swung away to the inside or outside by 90° (fig. 9), or be swung away and removed by pulling upward (fig. 10).

INFORMATION During reassembly, ensure that the footrest locking device has engaged.

The standard equipment includes removable quick-release axles on the rear wheels which facilitate transportation of your wheelchair.

To remove the rear wheels, press the button on the quick-release axle as described in Chapter 3. The rear wheels can then be removed. (fig. 1).

CAUTION **Accident risk.** During assembly, confirm that all bolts engage in the corresponding bore holes and that the quick-release rear wheel axles and casters securely lock into their axle housings. The rear wheels and casters must not be removable without the button being depressed! Check the function of the brakes.

4.1 Transfer

Transfer into the wheelchair is described below. The Motus is designed for persons who normally move around in the wheelchair independently. Therefore, the independent transfer is described. Should you require an attendant, the steps described here are carried out correspondingly but with support.

Maneuver the wheelchairs so that they are positioned at an angle of approx. 45° degrees to each other (fig. 11). Ensure that the wheel locks do not get caught in each other but remain freely movable. First activate the wheel locks of the wheelchair that is currently in use (fig.12).

Then activate the wheel locks of the wheelchair you want to transfer into (fig. 13). Flip up the footplates of the Motus (fig. 14). Move your feet from the footrest to the floor and move ahead in the wheelchair if required (fig. 15).

Transfer into the second wheelchair. The method illustrated is one possible way to support yourself (fig. 16). Find the most suitable path for you by practicing with an attendant.

Push the footplates of the Motus down again (fig. 14). Position your feet on the footplates of the wheelchair (fig. 17). Release the wheel locks (fig. 18). Now you can use the wheelchair.

5 Options

The wheelchair is designed as a modular assembly system. This means that certain accessories can be adapted to your wheelchair.

We recommend using a seat cushion at all times. The seat upholstery of the Motus has a loop texture so that seat cushions can be attached with hook and loop connection. All seat cushions listed on the order form have the required hook texture.

Furthermore, we will present you a selection of our variants and accessories which improve the use and comfort of your wheelchair.

5.1 Anti-tipper (Fig. 19/20)

Prevents the wheelchair from tipping too far backward. The anti-tipper can be adjusted in length without need of tools depending on seat height and seat angle (fig. 20). Strongly recommended for untrained wheelchair users as well as when extreme adjustments have been made to the rear wheel (for assembly instructions see Chapter 5.11).

5.2 Transport Wheels (Fig. 21)

With the rear wheels removed, the wheelchair can still be used as a rolling base; this makes it possible to manoeuvre through narrow passages (e.g. narrow bathroom doorways or airplane aisles).

⚠ CAUTION Accident risk. Please note that once the rear wheels are removed, the wheel locks can no longer be activated!

5.3 Crutch Holder (Fig. 22)

Allows crutches to be attached to the wheelchair.

5.4 Tip-assist (Fig. 23)

Makes it easier for an attendant to tip up the wheelchair, e.g. to cross a step. The tip-assist must be mounted so that the distance to the ground is at least **4 cm**.

5.5 Plug-on Side Panels (Fig. 24)

These side panels can be removed for getting into or out of the wheelchair. To remove the side panels, pull the out of the side panel attachment. Additionally, the armrest can be adjusted in height.

5.6 Desk Side Panels (Fig. 25)

The side panels can be swung away to the rear for transfer. For this purpose, press the half-moon push buttons and swing the side panel away to the rear. Remove the side panel in vertical position by pulling it upward. Moreover, the armrest is adjustable in height.

5.7 Plug-on Wheel Lock Lever Extension (Fig. 26)

The extension facilitates the wheel lock lever control to the users with limited hand function.


5.8 Rear Wheel with Drum Brake (Fig. 27)

The drum brakes allow the attendant to comfortably and safely brake the wheelchair. The rear wheels remain removable via the quick-release axle system.

5.9 Tray (Fig. 28)

Serves as a supporting surface during meals, when working or when playing. The clear material allows visibility of the legs and correction of the sitting posture.

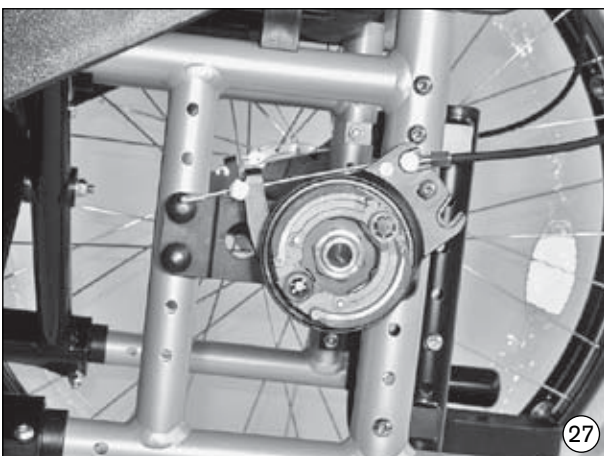5.10 Elevating Footrest (Fig. 29)

Allows positioning of the leg at different angles. Upon manipulation of the release lever which is integrated into the swivel segment the footrest angle can be adjusted in small increments.

5.11 Push Handles, Height-adjustable (Fig. 30)

Can be adjusted vertically to a more comfortable height for the attendant.


5.12 Wheelbase Extension (Fig. 31)

This option makes the wheelchair particularly stable against tipping and at the same time increases the turning radius.

INFORMATION For transfemoral amputees the long wheelbase setting is indispensable.

5.13 Angle Adjustable Back (Fig. 32)

The seat design allows the backrest angle adjustment by 30°. The unlocking of the adjustable back is activated via the pull cable.

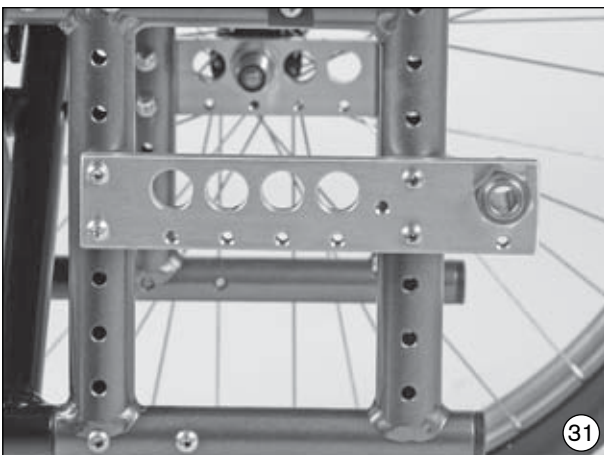5.14 Spoke Protector (not illustrated)

Protect the fingers from getting caught in the wheel spokes. Also improve the appearance of the wheelchair.

6 Adjustment / Assembly Instructions

INFORMATION

The following tools are required for repair and maintenance:

- Allen wrenches, sizes 4 and 5 mm
- Open-end wrenches, sizes: 10 mm, 11 mm, 13 mm, and 24 mm
- Set of screwdrivers
- Torque wrench
- Tire levers
- Tire pump

All wheelchair users have expectations regarding wheelchair type, accessories, seating and handling. Your Motus adaptive wheelchair offers many options that permit the wheelchair to be adapted according to your requirements. We recommend consulting with your dealer or therapist to determine which wheelchair settings will work best for you.

6.1 Adjusting Seat Height and Seat Inclination (Fig. 33)

The higher the rear wheel is attached to the frame, the more the seat inclines downward. This has two effects:

- First, the wheelchair has a greater tendency to tilt backwards.
- Secondly, the user sits deeper in the wheelchair and thus has more stability.


A careful and synchronized height adjustment of both the rear wheels and the casters will allow to adapt the seat height to the user's individual requirements.

The Motus offers the possibility to change vertical position of the rear wheel by shifting the rear wheel adapter along the frame. Loosen all four screw connections (Fig. 33, see arrow) and remove the screws. Shift the adapter to the desired position and firmly retighten the screws.

To achieve additional rear wheel positions, the rear wheel adapter can be turned by 180°.

INFORMATION When the rear wheel position is changed, the angle of the caster axle in relation to the ground also changes. Ensure that this angle is always approx. 90° wide (see Chapter 6.3). The wheel lock must be re-adjusted as well (see Chapter 6.4).

INFORMATION Make sure to firmly retighten all screws and nuts after making adjustments. The proper torque for the screw connection of the rear wheel adapter is **8 Nm**. The proper torque for the screw connection of the fitting is **50 Nm**.

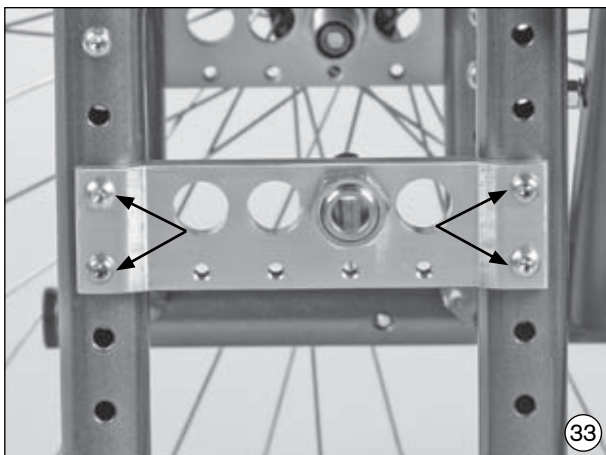6.2 Changing the Wheelbase (Fig. 34)

Shifting the rear wheel backwards lengthens the wheelbase and thus provides for greater stability of the wheelchair. Shifting the rear wheel forwards (only for advanced wheelchair drivers!) relieves the load on the casters. This increases the maneuverability of the wheelchair. The caster wheels can then also be lifted more easily to clear obstacles such as curbs or steps.

Loosen the screw connections from the fitting (fig. 34, see arrow) and remove the fitting. Shift the fitting to the desired position and firmly retighten the screws.

The distance between rear wheel and side panel can be continuously adjusted by adjusting the counter nut on the outer side of the rear wheel adapter.


When changing the adjustments originally made in the factory, two other component groups of the wheelchair must be readjusted, refer to Chapters 6.3 and 6.4.

⚠ CAUTION Risk of tipping. Please note: If the rear wheels are in a more forward position and the user's body is not appropriately positioned, the user may tip backward - even on level ground!

⚠ CAUTION Risk of tipping. For untrained wheelchair users and when there are extreme adjustments to the rear wheel, an anti-tipper is strongly recommended!

⚠ CAUTION Risk of tipping. For transfemoral amputees it is essential to shift the rear wheels backwards!

ℹ INFORMATION Make sure to firmly retighten all screws and nuts after making adjustments. The proper torque for the screw connection of the rear wheel adapter is **8 Nm**. The proper torque for the screw connection of the fitting is **50 Nm**.

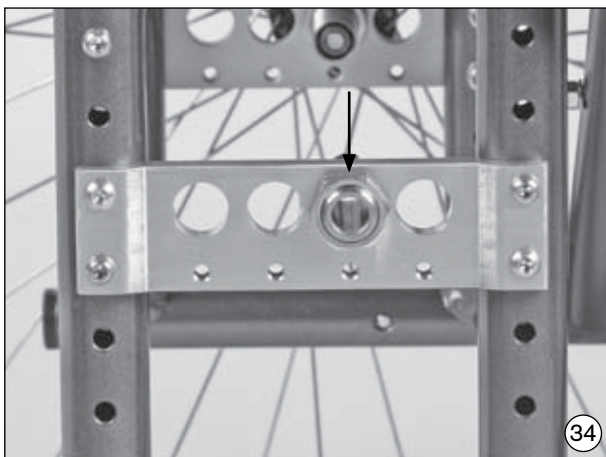6.3 Setting the Wheel Camber and Toe-in (Fig. 35-37)

The rear wheel adapter of the Motus allows for different camber angle settings of the rear wheels (0° and 2.5°). The greater the negative camber angle of the rear wheels, the more manoeuvrable the wheelchair becomes.

The rear wheel adapter itself has a wheel camber of 0° or 2.5°, depending on the arrangement of the wedge-shaped camber washers set at the factory. The wheel camber can be changed by turning round the mounted camber washers. Proceed as follows:

First disassemble the nuts on the rear wheel fitting.




Setting the rear wheel camber to 2.5°: Position the camber washer on the outside of the rear wheel adapter with the thinner side pointing upwards and into the desired engaging position (Fig. 35, item A). Place the second camber washer with the thicker side pointing upwards on the inside of the rear wheel adapter and into the corresponding engaging position (Fig. 36, item A).

Setting the rear wheel camber to 0°: Remove the camber washer from the outside of the rear wheel adapter and place it against the washer on the inside of the adapter such that the thinner side of the washers is next to the thicker side of the other washer (no Fig.).

Further proceeding: Once you have applied the retaining washer from the inside (Fig. 36, item B), slightly tighten the nuts on the rear wheel fitting. Check the space between the rear wheel and the side panel. You can adjust the space by rotating the fitting.

Toe-in correction: For toe-in correction, the nuts on the rear wheel fitting remain loosened. Then also loosen the oval head screw on the set plate (Fig. 35, item B). Now apply the level (provided) to the measuring point of the set plate and move the plate until it is level (Fig. 37). Then firmly retighten all screw connections.

INFORMATION Make sure to firmly retighten all screws and nuts after making adjustments. Observe the values indicated in section 6.1.

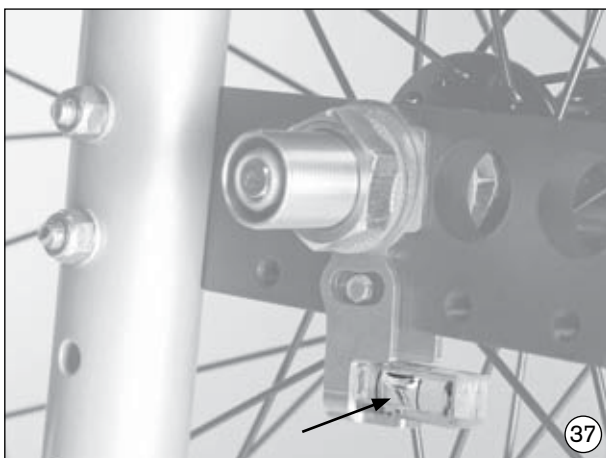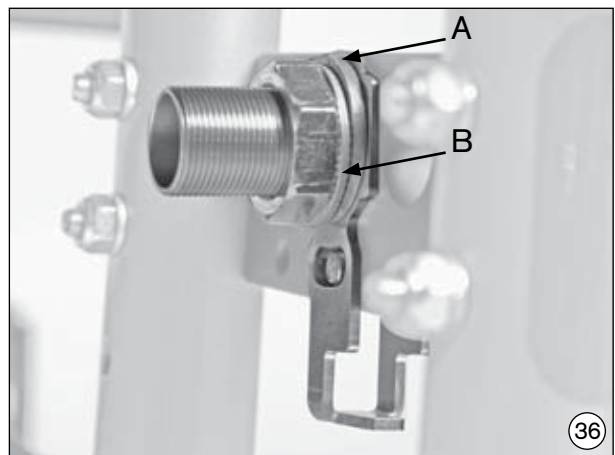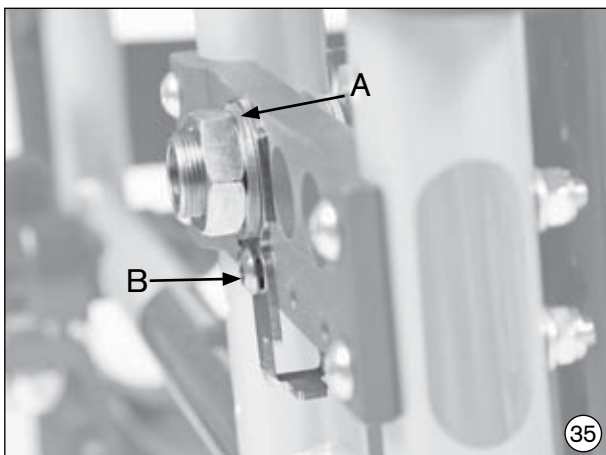6.4 Caster Trailing Adjustment (Fig. 38-41)

After the rear wheels have been mounted in the appropriate position, the caster wheel journal angle must be readjusted. The threaded axle should be as horizontal as possible in relation to the ground to ensure optimal driving behavior of the wheelchair. The caster adapter allows for continuous angle adjustment.

Remove the cap from caster head bearing with a screwdriver (fig. 38). Then remove the protective cap from the screw connection on the frame inner side (fig. 39).

Loosen the rear screw connection as shown in fig. 40.

Open the screw connection of the eccentric. Place the bubble level on the caster wheel journal and adjust the eccentric with a screwdriver (fig. 41). Retighten all screw connection and put on the cover caps.

INFORMATION The proper torque for tightening the screws is **8 Nm**.

6.5 Adjusting the Wheel Lock (Fig. 42)

To change the position of the rear wheels, first loosen the attachment screws of the clamping adapters of the wheel lock and push the wheel lock forward. Then reposition the rear wheels. The wheel lock should then be remounted such that the distance between the tires and wheel lock bolts is max. **5 mm** when the wheel lock is not activated (fig. 42; subject to technical changes).

⚠ CAUTION **Accident risk.** The wheel locks are only effective with sufficient air pressure and proper adjustment. The proper air pressure is indicated on the side wall of the tire. It should be at least **7 bar** (350 kPa). To ensure sufficient braking effect, use only original rear wheels with an approved radial excursion of maximally **± 1 mm**.

INFORMATION Make sure to firmly retighten all screws and nuts after making adjustments.

INFORMATION Tighten the wheel lock attachment with a torque of **10 Nm**.

6.6 Adjusting the Brake Force - Drum Brake (Fig. 43)

In order to achieve optimal braking effect, the brake force can be adjusted by means of the adjustment screw (fig. 43, item 1). The brake force is increased by turning the screw counter-clockwise and decreased by turning the screw clockwise.

Turn the adjustment screw counter-clockwise until frictional noise occurs at the turning wheel. Then tighten the adjustment screw until the frictional noise disappears. The wheel should run freely.

After finishing the adjustment, fasten the adjustment screw by tightening the counter nut.



INFORMATION Ensure that the brake force of both rear wheels is adjusted equally.

INFORMATION Make sure that the drum brake engages with sufficient force only when the manual brake lever is set to the second ratchet position.

6.7 Angle Adjustment of the Footplate (Fig. 44)

Prior to adjusting the footplate angle, loosen the screw connections as indicated in the figure (fig. 44, item A). Subsequently tilt the footplate until the desired angle is achieved.

INFORMATION Make sure to re-tighten the loosened screw connections. The proper torque is **6 Nm**.

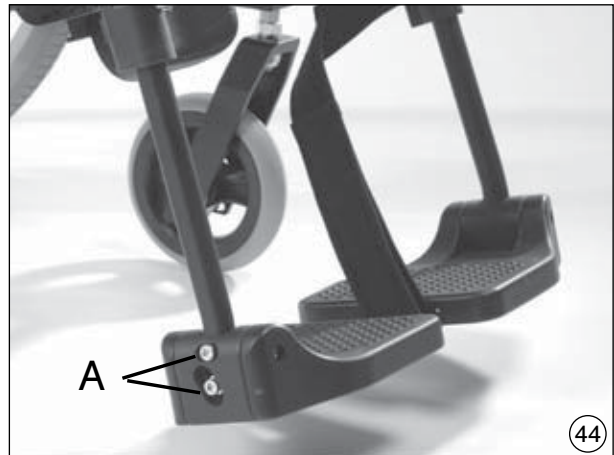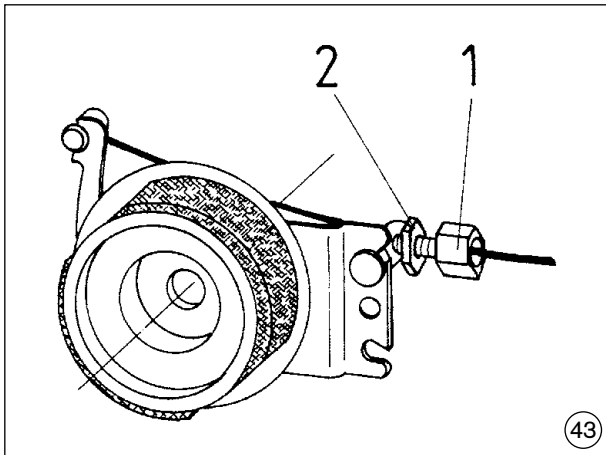6.8 Adjusting the Lower Leg Length (Fig. 45/46)

By loosening the set screw, the footplate can be adapted to your lower leg length and to the thickness of your seat cushion (standard footrest = fig. 45; elevating footrest = fig. 46).

Make sure when adjusting the footrest height, that the footrest bar is slid into the swivel segment at least **40 mm**. A marking on the footrest bar indicates how far the bar should be slid in.

INFORMATION After having made adjustments, retighten the set screw to a torque of **6 Nm**.



6.9 Adjusting the Angle of the Elevating Footrest (Fig. 47)

To adjust the footrest angle, turn the release lever (fig. 47, see arrow) up to the stop. While propping up the footrest bar, set the desired angle. Then carefully let the release lever turn back. At the next free position, the footrest will snap into place.

6.10 Adjusting the Height of the Armrest for Desk Side Panel (Fig. 48)

The armrest can be adjusted in height by bringing it into various positions without use of any tools. For this purpose, pull the button in the round opening in the side panel upward (fig. 48, see arrow) and shift the armrest into the desired position. The armrest automatically locks in place after releasing the button.

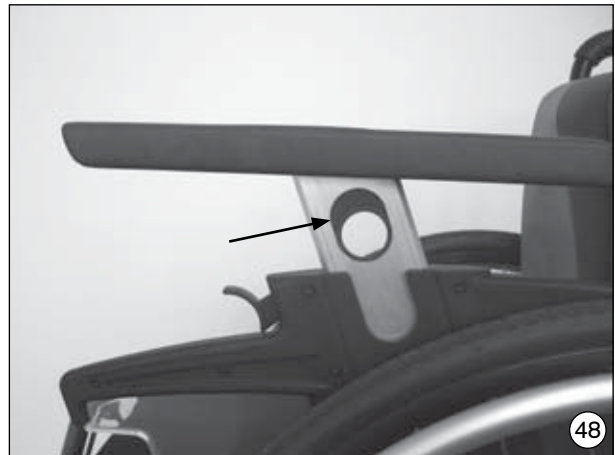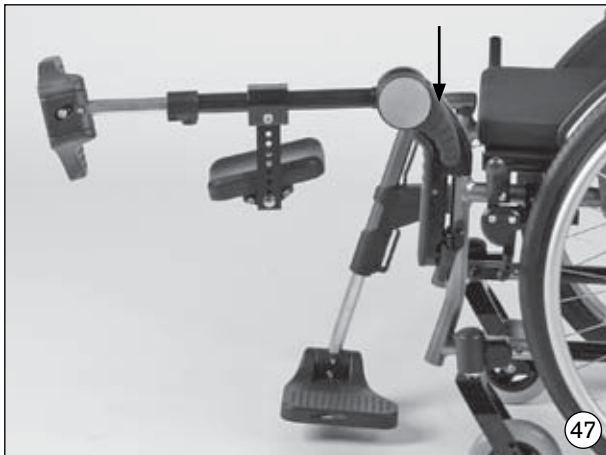6.11 Adjusting the Back Height (Fig. 49)

The back height is adjustable by **10 cm**. The back tube which is integrated into the posterior frame allows for the adjustment of the back height in increments of 2.5 cm.

For this purpose, remove the screws from the rear wheel adapter (fig. 49, see arrow) and set the desired back height.

INFORMATION Make sure to firmly retighten all screws and nuts after making adjustments. The proper torque is **8 Nm**.

6.12 Installing the Anti-tipper (Fig. 50–52)

1. Mount the receiver tube on the frame in the position as shown in the figure (fig. 50).
2. Screw the anti-tipper in the utmost position of the receiver tube as shown in the figure (fig. 51).
3. Adjust the anti-tipper using the longitudinal adjustment (fig. 52).
4. By changing the screw connection in the receiver tube the anti-tipper can be horizontally adjusted. The anti-tipper wheel must at least completely extend from the rear wheel to the back, and the distance between the wheel and ground should not exceed **5 cm** maximum. Find the appropriate position with the assistance of a helper!

6.13 Adapting the Back Upholstery (Fig. 53)

You can adjust the back upholstery in segments according to your needs.


To do this, remove the back upholstery padding and undo the hook-and-loop fasteners of the straps (fig. 53).

Reposition the straps as desired, attach them again with the hook-and-loop fasteners and pull on the back upholstery padding.

6.14 Changing the Back Angle (Fig. 54)

It may be useful to customize the preset angle between the seat and the backrest for therapeutic reasons. The angle can be adjusted between 90° and 120°.

Pull the release cable of the back release (fig. 54, arrow B) until the angle adjustment is released. Bring the backrest into the desired position and release the release cable. Ensure that the locking bolts are safely locked on both sides.

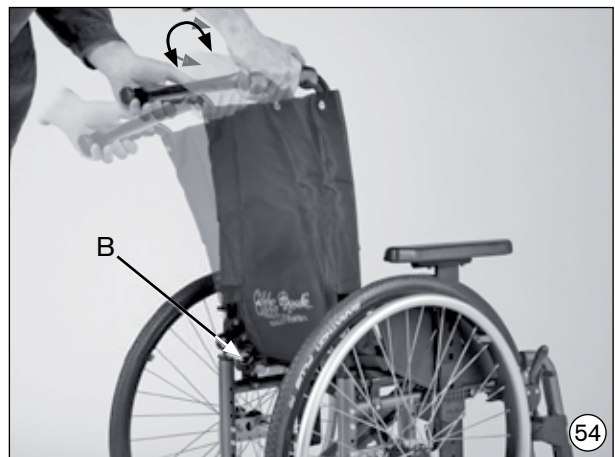6.15 Push Handles, Height-adjustable (Fig. 55)

To adjust the suitable height of the push handles, use the adjustable clamping lever to loosen the clamping screws.

INFORMATION Be sure to snugly re-tighten the clamping levers after all adjustments!

6.16 Attaching the Safety Belt (Fig. 56/57)

Attach the safety belt with the screws provided to the uppermost bore hole of the posterior frame part.

7 Maintenance, Cleaning, and Disinfection

7.1 Maintenance

Your wheelchair has been provided with the CE sign. The manufacturer herewith guarantees that this medical product as a whole conforms to the requirements of the European Directive for Medical Products 93/42/EEC.

The proper function of the wheelchair, especially of the brakes, should be checked before every use. Safety nuts should be used only once. If they have been loosened several times, they must be replaced.

Should any problems or wear be encountered, please contact your authorized specialist dealer immediately in order to have them corrected. This is particularly important if you notice that the driving characteristics change or the wheelchair becomes less stable. **We also recommend that you have the wheelchair serviced every 12 months by your Otto Bock authorized dealer.** The following table lists functions which the user must check at the indicated intervals.

Check	before each use	monthly	quarterly
Function test of the wheel locks and drum brake	X		
Sagging of the seat or back upholstery		X	
Footplate stability		X	
Sight-check of wearing parts (e.g. tires, bearings)		X	
Dirt on bearings		X	
Push ring damage		X	
Air pressure (see indication on tire)		X	
Wear on folding mechanism		X	
Spoke tension of the rear wheels			X
Screw connections			X

To ensure smooth operation at all times, the following precautions and cleanings should be carried out. Users may carry these out by themselves.

- Hair and dirt particles generally accumulate between the caster wheel and fork. This can restrict the caster wheels from rotating smoothly. Remove the caster and thoroughly clean the fork and caster using a mild detergent.
- The rear wheels feature a quick-release system. To keep this system operational, ensure that no dirt adheres to the quick-release axle or axle housing. The quick-release axle should also be lightly lubricated regularly with resin-free sewing machine oil.
- If your wheelchair gets wet, towel dry it as soon as possible.
- Do not use your wheelchair in salt water and keep sand or other particles from damaging the wheel bearings.

- Screw connections should be checked frequently, in particular when beginning to use the wheelchair and after any adjustment. If a screw connection becomes loose repeatedly, consult your dealer.

7.2 Cleaning and Disinfection

Clean the padding and upholstery with warm water and a washing-up liquid. Remove spots with a sponge or a soft brush. Rinse with clear water and let the components dry.

INFORMATION Do not use any aggressive cleansing agents, solvents, or hard brushes etc.

INFORMATION No wet-cleaning. No machine-wash.

For disinfection, water-based agents (e.g. Sagrothan original concentrate) should be used. The manufacturer's instructions for use must be observed.

INFORMATION Before disinfecting, the pads and handles must be cleaned.

Plastic parts, frame parts as well as chassis and wheels can be cleaned with a mild detergent. Dry thoroughly.

8 Changing Tires (Fig. 58/59)

Changing flat tires requires only the necessary tools and minimal technical skills. Users may change tires themselves if they wish. We recommend carrying an emergency repair set and an air pump at all times (if using pneumatic tires).

Suitable tire pumps are listed on the order sheet and are supplied with the product. Alternatively, users can keep a spray can of hardening foam (available at bicycle shops) to fill a tire in case of punctures.

- In case of a flat tire, carefully remove the tire from the rim using the appropriate tools.
- Take care not to damage the rim or the inner tube.
- Repair the inner tube according to the directions in the repair kit, or, replace the damaged inner tube with a new one.
- Before re-mounting the tire, inspect the inner rim surface and the inside of the tire for any object that may have caused the puncture.
- Be sure that the rim bands are not damaged. They protect the tube from being damaged by the ends of the spokes.

Mounting (Fig. 58)

- Place the protective rim band over the valve. Then push the valve through the rim. Screw on the valve nut, if there is one. The rim band can then be easily applied.
- Ensure that all spoke ends are covered.

Tires (Fig. 59)

- Beginning behind the valve, place the lower side of the tire over the edge of the wheel rim. Pump the tube with a small amount of air until it assumes its round shape. Then insert it into the tire.
- Check the inner tube for folds. If folds are present, release some air. Now you can easily install the upper side of the tire starting across from the valve, using both hands and working towards the valve.

Inflating


- Check around the tire on both sides to ensure that the tube is not clamped between the edge of the tire and the rim.
- Push the valve back slightly and then pull it out again such that the tire is well-seated in the valve area.
- Inflate the tire until it can still be pressed easily with your thumb. If the control line on both sides of the tire is equally distanced to the rim, the tire sits centrally. If the tire does not sit centrally, let out some air and realign the tire. Then inflate the tire to the pressure indicated on the side wall and tighten the dust cap.

9 Technical Data (in cm)

Seat width:	35.5 – 55.5
Seat depth:	36.0 – 54.0
Seat angle:	adjustable from -5° to 15°
Seat height:	42.5 – 49.0 caster fork short
	46.0 – 53.0 caster fork long
Back angle:	adjustable up to -30°
Lower leg length:	34.0 – 54.0
Armrest height:	21.0 – 30.0
Anterior support area of the armrest:	30.0
Overall width:	52.0 – 72.0 with standard rear wheels*
	55.5 – 75.5 with rear wheels with drum brake*
Overall length:	80.5 – 108.5 with 22" rear wheel
	83.0 – 111.0 with 24" rear wheel
Horizontal position of the rear wheel axle:	5.0
Overall width, folded	32.0 – 35.5
Overall length, folded:	80.5 – 111.0
Overall height, folded:	73.0 – 109.0
Weight:	approx. 13 kg (29 lbs)
Static stability in forward direction:	10°
Static stability in sideward direction:	10°
Static stability in rearward direction:	0 - 13°, depending on adjustment. Observe safety instructions!
Weight of the heaviest part:	from 8 kg (17.6 lbs)
Load capacity:	125 / 140 kg (275 / 308 lbs)


* Applicable to push ring attachment, narrow (for push ring attachment, wide: + 20 mm) and a 0° camber of the rear wheels.

Anterior seat height (in cm)

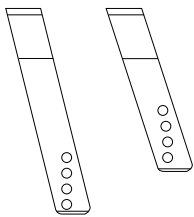	Item	Caster diameter				
		4"	5"	140 mm	6"	7"
Caster fork, short (481F00-SS072)	4					
	3	42.5				
	2	43.5	45.0	46.0	46.0	
	1	45.5	46.0	47.0	48.0	49.0
Caster fork, long (481F00-SS070)	4		46.0	47.0	47.5	49.0
	3	46.0	47.5	48.0	49.0	50.0
	2	47.5	49.0	50.0	50.0	51.0
	1	48.5	50.0	51.5	51.0	53.0

Posterior seat height (in cm)

Item	Rear wheel size	
	22"	24"
1	36.0	38.5
2	37.0	40.0
3	38.5	41.0
4	40.0	42.5
5	41.0	43.5
6	42.0	44.5
7	43.5	46.0
8	44.5	47.5
9	46.0	48.5
10	47.0	49.5
11	48.5	51.0
12	50.0	52.0

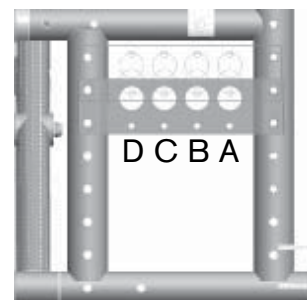INFORMATION The values indicated are measurements which have been theoretically determined. Not all adjustment possibilities can be used for all wheelchair models. Further, the adjustment combinations are limited due to the sturdy frame geometry. We reserve the right to technical alterations and deviations by 5%.

10 Information on Re-use

CAUTION

Risk due to improper use. Seat padding which come into direct contact with the skin may cause functional or hygienic risks for another patient using one and the same product. They should be replaced in case of re-use.

The product is suitable for re-use.

Similar to second-hand machines or cars, products that are being re-used are subject to increased strain. Features and functions must not change in a way that could endanger patients or other persons within the product's life cycle.

Based on market observations and the current state of technology, the manufacturer has calculated that the product can be used for a period of **4 years**, provided that it is used properly and that the service and maintenance instructions are observed. Periods during which the wheelchair is stored at the dealer or with the third party payer are not included in this period. It should be clearly pointed out, however, that the product is reliable far beyond this defined period of time, provided that it is cared for and maintained appropriately.

In cases of re-use, the corresponding product must first be thoroughly cleaned and disinfected. Subsequently, the product must be examined by an authorized specialist to check the condition, wear, and possible damages. Any worn and damaged components as well as components which do not fit or are not suitable for the user must be replaced.

The service manual includes a service schedule, detailed information, and a list of the required tools.

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Otto Bock Healthcare Products GmbH

Kaiserstraße 39 · A-1070 Wien
Tel. +43 1 5269548 · Fax +43 1 5267985
e-mail: vertrieb.austria@ottobock.com
www.ottobock.at

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Otto Bock HealthCare LP

Two Carlson Parkway North, Suite 100
U.S.A.-Minneapolis, Minnesota 55447
Phone +1 800 328 4058 · Fax +1 800 962 2549
e-mail: usa.customerservice@ottobock.com
www.ottobockus.com

Asia/Pacific

Otto Bock Australia Pty. Ltd.

Suite 1.01, Century Corporate Centre
62 Norwest Boulevard · Norwest Business Park
AUS-Baulkham Hills NSW 2153
Tel. +61 2 88182800 · Fax +61 2 88182898
e-mail: healthcare@ottobock.com.au
www.ottobock.com.au

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Otto Bock Mobility Solutions GmbH has been certified by the German Society for the Certification of Quality Assurance Systems (DQS) in accordance with DIN EN ISO 9001 standard, reg. no. 779 (management system)

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Otto Bock Korea HealthCare Inc.
Otto Bock South East Asia Co. Ltd.
Otto Bock HealthCare GmbH
Otto Bock Manufacturing Königsee GmbH